happylegs®


Instrucciones
User manual
Mode d'emploi
Bedienungsanleitung
Manuale d'uso
Manual do usuário

www.happylegs.es

happylegs

la máquina de andar sentado

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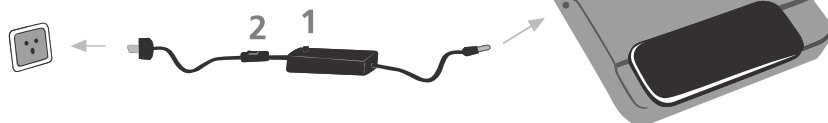
the seated walking machine

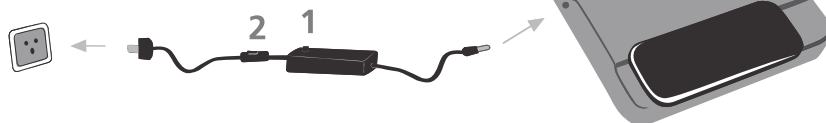
Getting Started:

The box contains **1 Happylegs** and **1 Power Transformer***.

The powered exerciser Happylegs has an energy-efficient motor of low power consumption. The equipment connects to the mains through 3-speed transformer. To use Happylegs simply connect the transformer to the connector at the front of the Happylegs machine and then plug the transformer socket into mains electricity. To put on and off the machine you have to use the switch in the cable (2). You can select the required speed through the speed selector (1).

You can lean the Happylegs unscrewing the 2 height-adjustable legs at the lower side (3). Do not force the motion because Happylegs automatically moves your legs without any effort.

**The Happylegs Premium model includes a remote control kit. See instructions inside.*

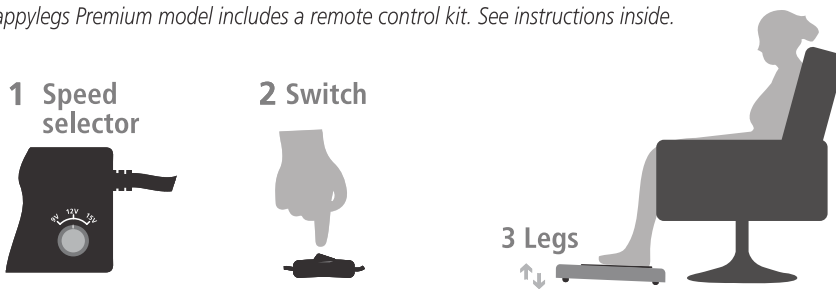Select the required operational speed:

Turn on the transformer (2) and select the most suitable speed (1).

- **9 V Slowest:** This speed is the more gentle an ideal for people who may be recovering from injuries or have stiff joints and for those who suffer fatigue. It can be used during hours.
- **12 V Normal:** Recommended for prolonged use, around 6-8 hours at home or at work.
- **15 V Fast:** The use of this speed depends on the age and health of the user. If the user is elderly or suffers from fatigue the recommended use is for periods of 2 to 3 hours daily.

The rest of the users can use it between 6- 8 hours.

Cleaning and Maintenance:

Never stand on top of the machine, as full weight could damage the internal operational mechanism and it is out of the warranty of the product. Happylegs does not require regular maintenance and has been developed to be used continuously for long periods of time. In order to clean it, a dry cloth should be used when the equipment is unplugged and switched off.

Transportation and Storage:

As this product is electrically operated, avoid spilling any liquid on the Happylegs or storing in high humidity environments. Also avoid temperatures greater than 60° C. Happylegs should be transported within its packaging. In case of failure, please contact your official supplier. Not to be disassembled or repaired by unauthorized personnel.

Technical specifications:

Happylegs is powered by mains electricity via a transformer designed exclusively for Autoergon 2002. Do not use parts or accessories not supplied by the official Happylegs manufacturer or supplier.

- Connected to the mains power.
- Input Voltage: 110-240 V// 50-60Hz 1,5A 
- Output Voltage: 9V/ 12V/ 15V 3A 

B  .  Class II.    RoHS

Warranty Terms and CE Marking

Happylegs has been manufactured with the highest quality both in the materials used and quality control and design of its components to obtain the European Quality certification (CE Marking). The warranty terms covers all possible faults caused by manufacturing defects or defects in the materials used in the manufacturing of the machine. The duration of the warranty of the machine Happylegs is 5 years from the date of the purchase. The warranty of external components is 2 years. If any component fail to work you have to contact your official supplier in the country where you have purchased it.

The warranty does not cover you when:

- The machine is used incorrectly or in a different way from the operating instructions.
- Breakages have been caused by rough handling, knocking, dropping or crushing, destruction by fire or flooding.
- The power cable fails or breaks down for misuse or negligence.
- Dismantling the machine or handling it by other people that do not work for Autoergon 2002.


The holder of the warranty is in full possession of the minimum rights recognised by law (Law 23/2003 on Guarantees in the Sale of Consumer Goods, published in B.O.E. - State Official Bulletin - the 11th July 2003, that includes the 1999/44/CE Directive of the European Parliament and the Council, of 25th May 1999).

Disposal and recycling of electrical and electronic appliances:

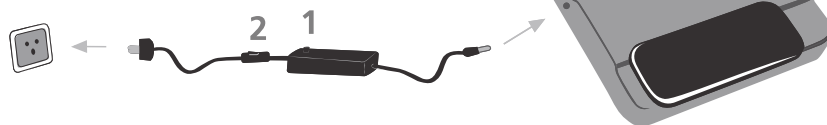1. If the symbol with a crossed container is displayed, this means it complies with Directive 2002/96/CE.
2. All electrical or electronic appliances should be disposed properly at designated located places only. The correct disposal and treatment of damaged goods contributes to avoid potential risks to the health and the environment.
3. To obtain more information about disposal of your electrical and/or electronic devices, contact your local authorities or the retail place of purchase.

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